


MERCK

TRUST OUR NMR expertise

One-Stop Shop for NMR Consumables

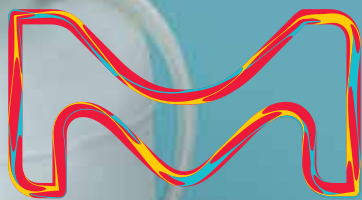

The life science
business of Merck
operates as
MilliporeSigma in
the U.S. and Canada.

Sigma-Aldrich[®]
Lab & Production Materials

Not just
good or great.
But gold.

As a leader in life sciences, Merck offers a catalog of more than 300,000 products. One of our core specialties is NMR, for which we produce premium quality solvents, standards, CRMs, and stable isotope reagents, as well as dedicated NMR tubes and accessories.

Our innovative certified reference materials (CRMs) are ideal for a wide range of new applications in quantitative NMR (qNMR), such as metabolomics, biomarker discovery, and physiological pathway research. All our CRMs are traceable to NIST SRM, and produced according to ISO/IEC 17025 and ISO Guide 34 criteria. This double accreditation represents the “gold standard” in reference materials – the highest quality possible.

For more information, please visit
[SigmaAldrich.com/NMR](https://www.sigmaaldrich.com/NMR)

TABLE OF CONTENTS:

High Quality NMR Solvents	3
TraceCERT® Certified Reference Materials for quantitative NMR	9
NMR Reference Standards for Calibration	12
Stable Isotope Labeled Reagents and Building Blocks	13
NMR Tubes and Accessories	14



High Quality NMR Solvents

High quality solvents are essential for the rigorous demands of NMR research and analysis. At Merck, we're passionate about quality. That's why we provide the widest range of NMR solvents with excellent chemical purity and the highest isotopic enrichment possible.

To make sure we always offer you the finest quality, we consistently review and improve our methods for solvent purification and water content reduction. And, all our NMR solvents undergo strict quality control testing during production and packaging.

For more specific requirements, we support you with the convenience of customized products and packaging. Our glass bottles, ampules, and septum bottles for NMR solvents are available in the broadest range of configurations and sizes.

Benefits

- Reliable results in NMR spectra
- Excellent chemical purity and highest isotopic enrichment
- Reliable deuteration degrees
- Determination of water content in Karl Fischer and NMR
- Safe, easy handling with septum bottles and glass ampules
- Flexible choice of packaging formats

Cat. No.	Product Description	Isotopic Purity and Internal Standards
Acetic acid-d		
151777-50G	Acetic acid-d	99 atom % D
151777-250G	Acetic acid-d	99 atom % D
151785-50G	Acetic acid-d ₄	≥99.5 atom % D
151785-10X1ML	Acetic acid-d ₄	≥99.5 atom % D
151785-25G	Acetic acid-d ₄	≥99.5 atom % D
151785-10G	Acetic acid-d ₄	≥99.5 atom % D
151785-5G	Acetic acid-d ₄	≥99.5 atom % D
416886-25G	Acetic acid-d ₄	≥99.5 atom % D, contains 0.03 % (v/v) TMS
233315-5G	Acetic acid-d ₄	≥99.9 atom % D
233315-10X0.5ML	Acetic acid-d ₄	≥99.9 atom % D
Acetone-d₆		
151793-100G	Acetone-d ₆	99.9 atom % D
151793-25G	Acetone-d ₆	99.9 atom % D
151793-10X1ML	Acetone-d ₆	99.9 atom % D
151793-10X0.75ML	Acetone-d ₆	99.9 atom % D
151793-10G	Acetone-d ₆	99.9 atom % D
151793-50G	Acetone-d ₆	99.9 atom % D
444863-1L	Acetone-d ₆	99.9 atom % D
444863-10ML-AMP	Acetone-d ₆	99.9 atom % D
444863-25ML	Acetone-d ₆	99.9 atom % D
444863-100ML	Acetone-d ₆	99.9 atom % D
444863-10ML-GL	Acetone-d ₆	99.9 atom % D
444863-10X0.6ML	Acetone-d ₆	99.9 atom % D
444863-10X0.5ML	Acetone-d ₆	99.9 atom % D
444863-50ML	Acetone-d ₆	99.9 atom % D
434531-10X1ML	Acetone-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
434531-25G	Acetone-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
434531-10G-GL	Acetone-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
434531-50G	Acetone-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
454133-50G	Acetone-d ₆	99.9 atom % D, contains 1 % (v/v) TMS
454133-25G	Acetone-d ₆	99.9 atom % D, contains 1 % (v/v) TMS

Cat. No.	Product Description	Isotopic Purity and Internal Standards
175862-10X0.75ML	Acetone-d ₆	99.96 atom % D
175862-5G	Acetone-d ₆	99.96 atom % D
175862-10X0.5ML	Acetone-d ₆	99.96 atom % D
Acetonitrile-d₃		
151807-50G	Acetonitrile-d ₃	≥99.8 atom % D
151807-5G	Acetonitrile-d ₃	≥99.8 atom % D
151807-10G-GL	Acetonitrile-d ₃	≥99.8 atom % D
151807-10X1ML	Acetonitrile-d ₃	≥99.8 atom % D
151807-10X0.6ML	Acetonitrile-d ₃	≥99.8 atom % D
151807-25G	Acetonitrile-d ₃	≥99.8 atom % D
151807-10G-AMP	Acetonitrile-d ₃	≥99.8 atom % D
151807-10X0.75ML	Acetonitrile-d ₃	≥99.8 atom % D
151807-100G	Acetonitrile-d ₃	≥99.8 atom % D
699543-10G	Acetonitrile-d ₃	≥99.8 atom % D, Special HOH
699543-50G	Acetonitrile-d ₃	≥99.8 atom % D, Special HOH
569550-10X1ML	Acetonitrile-d ₃	≥99.8 atom % D, anhydrous
366544-5G	Acetonitrile-d ₃	99.8 atom % D, contains 0.03 % (v/v) TMS
366544-10X0.6ML	Acetonitrile-d ₃	99.8 atom % D, contains 0.03 % (v/v) TMS
366544-50G	Acetonitrile-d ₃	99.8 atom % D, contains 0.03 % (v/v) TMS
366544-25G	Acetonitrile-d ₃	99.8 atom % D, contains 0.03 % (v/v) TMS
233331-50G	Acetonitrile-d ₃	99.8 atom % D, contains 1 % (v/v) TMS
233323-10X0.5ML	Acetonitrile-d ₃	99.96 atom % D
233323-10X0.75ML	Acetonitrile-d ₃	99.96 atom % D
233323-5G	Acetonitrile-d ₃	99.96 atom % D
233323-25G	Acetonitrile-d ₃	99.96 atom % D
Benzene-d₆		
175978-10G	Benzene-d ₆	99 atom % D
175978-50G	Benzene-d ₆	99 atom % D
151815-10G-AMP	Benzene-d ₆	99.6 atom % D
151815-10G-GL	Benzene-d ₆	99.6 atom % D
151815-100G	Benzene-d ₆	99.6 atom % D
151815-10X1ML	Benzene-d ₆	99.6 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
151815-10X0.75ML	Benzene-d ₆	99.6 atom % D
151815-25G	Benzene-d ₆	99.6 atom % D
151815-50G	Benzene-d ₆	99.6 atom % D
570680-50G	Benzene-d ₆	anhydrous, ≥99.6 atom % D
364940-50G	Benzene-d ₆	99.6 atom % D, contains 0.03 % (v/v) TMS
364940-10G	Benzene-d ₆	99.6 atom % D, contains 0.03 % (v/v) TMS
364940-25G	Benzene-d ₆	99.6 atom % D, contains 0.03 % (v/v) TMS
561509-10X.75ML	Benzene-d ₆	99.96 atom % D, contains 0.03 % (v/v) TMS
175870-25G	Benzene-d ₆	99.96 atom % D
175870-10X0.5ML	Benzene-d ₆	99.96 atom % D
175870-10X1ML	Benzene-d ₆	99.96 atom % D
Bromobenzene-d₅		
175730-5G	Bromobenzene-d ₅	99.5 atom % D
175730-25G	Bromobenzene-d ₅	99.5 atom % D
tert-Butanol-d₁₀		
175889-1G	tert-Butanol-d ₁₀	99 atom % D
175889-5G	tert-Butanol-d ₁₀	99 atom % D
tert-Butyl methyl-d₃ ether		
434132-1G	tert-Butyl methyl-d ₃ ether	99 atom % D
Chlorobenzene-d₅		
176605-5G	Chlorobenzene-d ₅	99 atom % D
176605-1G	Chlorobenzene-d ₅	99 atom % D
Chloroform-d		
151823-10X0.75ML	Chloroform-d	99.8 atom % D
151823-250G	Chloroform-d	99.8 atom % D
151823-1KG	Chloroform-d	99.8 atom % D
151823-100G	Chloroform-d	99.8 atom % D
151823-50G	Chloroform-d	99.8 atom % D
151823-10X0.6ML	Chloroform-d	99.8 atom % D
151823-10X1ML	Chloroform-d	99.8 atom % D
570699-50G	Chloroform-d	≥99.8 atom % D, anhydrous
225789-10X0.6ML	Chloroform-d	99.8 atom % D, contains 0.03 % (v/v) TMS
225789-100G	Chloroform-d	99.8 atom % D, contains 0.03 % (v/v) TMS
225789-500G	Chloroform-d	99.8 atom % D, contains 0.03 % (v/v) TMS
225789-10X1ML	Chloroform-d	99.8 atom % D, contains 0.03 % (v/v) TMS
225789-50G	Chloroform-d	99.8 atom % D, contains 0.03 % (v/v) TMS
612200-100G	Chloroform-d	99.8 atom % D, contains 0.05 % (v/v) TMS
434876-500G	Chloroform-d	99.8 atom % D, contains 0.1 % (v/v) TMS
434876-150G	Chloroform-d	99.8 atom % D, contains 0.1 % (v/v) TMS
434876-100G	Chloroform-d	99.8 atom % D, contains 0.1 % (v/v) TMS
612200-100G	Chloroform-d	99.8 atom % D, contains 0.05 % (v/v) TMS

Cat. No.	Product Description	Isotopic Purity and Internal Standards
151831-50G	Chloroform-d	99.8 atom % D, contains 1 % (v/v) TMS
151831-100G	Chloroform-d	99.8 atom % D, contains 1 % (v/v) TMS
151831-250G	Chloroform-d	99.8 atom % D, contains 1 % (v/v) TMS
416754-250G	Chloroform-d	≥99.8 atom % D, contains 0.5 wt. % silver foil as stabilizer
416754-100G	Chloroform-d	≥99.8 atom % D, contains 0.5 wt. % silver foil as stabilizer
530735-250G	Chloroform-d	≥99.8 atom % D, contains 0.5 wt. % silver foil as stabilizer, 0.03 % (v/v) TMS
530735-100G	Chloroform-d	≥99.8 atom % D, contains 0.5 wt. % silver foil as stabilizer, 0.03 % (v/v) TMS
151858-10X0.75ML	Chloroform-d	99.96 atom % D
151858-50G	Chloroform-d	99.96 atom % D
151858-10X0.5ML	Chloroform-d	99.96 atom % D
151858-10X0.25ML	Chloroform-d	99.96 atom % D
151858-10X1ML	Chloroform-d	99.96 atom % D
151858-10G	Chloroform-d	99.96 atom % D
431915-10ML	Chloroform-d	99.96 atom % D, contains 0.5 wt. % silver wire as stabilizer
431915-50ML	Chloroform-d	99.96 atom % D, contains 0.5 wt. % silver wire as stabilizer
494275-10X0.75ML	Chloroform-d	99.96 atom % D, contains 0.03 % (v/v) TMS
494275-50G	Chloroform-d	99.96 atom % D, contains 0.03 % (v/v) TMS
494275-10G	Chloroform-d	99.96 atom % D, contains 0.03 % (v/v) TMS
Cyclohexane-d₁₂		
151866-10X1ML	Cyclohexane-d ₁₂	≥99.6 atom % D
151866-5G	Cyclohexane-d ₁₂	≥99.6 atom % D
151866-10G	Cyclohexane-d ₁₂	≥99.6 atom % D
151866-1G	Cyclohexane-d ₁₂	≥99.6 atom % D
Decahydronaphthalene-d₁₈		
217131-1G	Decahydronaphthalene-d ₁₈	98 atom % D
Deuterium oxide		
435767-100G	Deuterium oxide	99 atom % D
435767-25G	Deuterium oxide	99 atom % D
435767-1KG	Deuterium oxide	99 atom % D
617385-1KG	Deuterium oxide	99.8 atom % D
617385-1.107KG	Deuterium oxide	99.8 atom % D
151882-1.107KG	Deuterium oxide	99.9 atom % D
151882-10X1ML	Deuterium oxide	99.9 atom % D
151882-500G	Deuterium oxide	99.9 atom % D
151882-10X0.6ML	Deuterium oxide	99.9 atom % D
151882-100G	Deuterium oxide	99.9 atom % D
151882-250G	Deuterium oxide	99.9 atom % D
151882-1KG	Deuterium oxide	99.9 atom % D
151882-25G	Deuterium oxide	99.9 atom % D
151882-10X0.75ML	Deuterium oxide	99.9 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
151882-4KG	Deuterium oxide	99.9 atom % D
151882-10G	Deuterium oxide	99.9 atom % D
613444-100G	Deuterium oxide	99.9 atom % D, glass distilled
293040-100G	Deuterium oxide	99.9 atom % D, contains 0.75 wt. % 3-(trimethylsilyl) propionic-2,2,3,3-d4 acid, sodium salt
293040-25G	Deuterium oxide	99.9 atom % D, contains 0.75 wt. % 3-(trimethylsilyl) propionic-2,2,3,3-d4 acid, sodium salt
343773-25G	Deuterium oxide	99.9 atom % D, contains 1 % (w/w) 3-(trimethylsilyl)-1-propanesulfonic acid, sodium salt (DSS)
343773-100G	Deuterium oxide	99.9 atom % D, contains 1 % (w/w) 3-(trimethylsilyl)-1-propanesulfonic acid, sodium salt (DSS)
450510-10X0.75ML	Deuterium oxide	99.9 atom % D, contains 0.05 wt. % 3-(trimethylsilyl) propionic-2,2,3,3-d4 acid, sodium salt
450510-100ML	Deuterium oxide	99.9 atom % D, contains 0.05 wt. % 3-(trimethylsilyl) propionic-2,2,3,3-d4 acid, sodium salt
450510-25ML	Deuterium oxide	99.9 atom % D, contains 0.05 wt. % 3-(trimethylsilyl) propionic-2,2,3,3-d4 acid, sodium salt

Cat. No.	Product Description	Isotopic Purity and Internal Standards
296163-10G	Dichloromethane-d ₂	≥99.5 atom % D, contains 0.03 % (v/v) TMS
296163-1G	Dichloromethane-d ₂	≥99.5 atom % D, contains 0.03 % (v/v) TMS
296163-25G	Dichloromethane-d ₂	≥99.5 atom % D, contains 0.03 % (v/v) TMS
444324-25G	Dichloromethane-d ₂	99.9 atom % D
444324-10X0.75ML	Dichloromethane-d ₂	99.9 atom % D
444324-10X1ML	Dichloromethane-d ₂	99.9 atom % D
444324-10X0.6ML	Dichloromethane-d ₂	99.9 atom % D
444324-10G	Dichloromethane-d ₂	99.9 atom % D
444324-5G	Dichloromethane-d ₂	99.9 atom % D
530506-25G	Dichloromethane-d ₂	99.9 atom % D, contains 0.1 % (v/v) TMS
530506-1G	Dichloromethane-d ₂	99.9 atom % D, contains 0.1 % (v/v) TMS
233366-10X0.5ML	Dichloromethane-d ₂	99.96 atom % D
233366-5G	Dichloromethane-d ₂	99.96 atom % D
N,N-Dimethylacetamide-d₉		
522414-5G	N,N-Dimethylacetamide-d ₉	99 atom % D
522414-1G	N,N-Dimethylacetamide-d ₉	99 atom % D
N,N-Dimethylformamide-d₇		
189979-5X0.5ML	N,N-Dimethylformamide-d ₇	≥99.5 atom % D
189979-1G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D
189979-10G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D

Safe, easy handling with septum bottles and glass ampules

Cat. No.	Product Description	Isotopic Purity and Internal Standards
151890-50G	Deuterium oxide	≥99.96 atom % D
151890-250G	Deuterium oxide	≥99.96 atom % D
151890-10G	Deuterium oxide	≥99.96 atom % D
151890-10X0.75ML	Deuterium oxide	≥99.96 atom % D
151890-10X1ML	Deuterium oxide	≥99.96 atom % D
151890-1KG	Deuterium oxide	≥99.96 atom % D
151890-125G	Deuterium oxide	≥99.96 atom % D
191701-10X0.5ML	Deuterium oxide	99.990 atom % D
191701-10X0.25ML	Deuterium oxide	99.990 atom % D
191701-10G	Deuterium oxide	99.990 atom % D
613398-50G	Deuterium oxide	extra, 99.994 atom % D
613398-10G	Deuterium oxide	extra, 99.994 atom % D
Dibromomethane-d₂		
259020-5G	Dibromomethane-d ₂	99 atom % D, contains copper as stabilizer
177865-10X1ML	Dichloromethane-d ₂	99.5 atom % D
177865-25G	Dichloromethane-d ₂	99.5 atom % D
177865-10G	Dichloromethane-d ₂	99.5 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
189979-5G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D
189979-10X0.75ML	N,N-Dimethylformamide-d ₇	≥99.5 atom % D
189979-10X1ML	N,N-Dimethylformamide-d ₇	≥99.5 atom % D
269905-1G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D, contains 1 % (v/v) TMS
269905-10G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D, contains 1 % (v/v) TMS
700428-10G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D, contains 0.03 % (v/v) TMS
700428-1G	N,N-Dimethylformamide-d ₇	≥99.5 atom % D, contains 0.03 % (v/v) TMS
Dimethyl sulfide-d₆		
416452-1G	Dimethyl sulfide-d ₆	99 atom % D
175943-100G	Dimethyl sulfoxide-d ₆	99.5 atom % D
175943-50G	Dimethyl sulfoxide-d ₆	99.5 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
175943-250G	Dimethyl sulfoxide-d ₆	99.5 atom % D
175943-10G	Dimethyl sulfoxide-d ₆	99.5 atom % D
151874-25G	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10G-GL	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-50G-SB	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10X0.75ML	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10G-AMP	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-100G	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-5G	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-40KG	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-600G	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-50G-GL	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10X0.6ML	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10G-SB	Dimethyl sulfoxide-d ₆	99.9 atom % D
151874-10X1ML	Dimethyl sulfoxide-d ₆	99.9 atom % D
716731-50ML	Dimethyl sulfoxide-d ₆	≥99.9 atom % D, "Special HOH"
716731-10ML	Dimethyl sulfoxide-d ₆	≥99.9 atom % D, "Special HOH"
716731-10X0.75ML	Dimethyl sulfoxide-d ₆	≥99.9 atom % D, "Special HOH"
570672-10X1ML	Dimethyl sulfoxide-d ₆	anhydrous, 99.9 atom % D
570672-50G	Dimethyl sulfoxide-d ₆	anhydrous, 99.9 atom % D
296147-10X1ML	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-10X0.75ML	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-10X0.5ML	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-10X0.6ML	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-10G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-50G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
296147-25G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 0.03 % (v/v) TMS
185965-50G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 1 % (v/v) TMS
185965-25G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 1 % (v/v) TMS
185965-10G	Dimethyl sulfoxide-d ₆	99.9 atom % D, contains 1 % (v/v) TMS
156914-10X0.75ML	Dimethyl sulfoxide-d ₆	99.96 atom % D
156914-5G	Dimethyl sulfoxide-d ₆	99.96 atom % D
156914-10X0.5ML	Dimethyl sulfoxide-d ₆	99.96 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
156914-10X1ML	Dimethyl sulfoxide-d ₆	99.96 atom % D
156914-1G	Dimethyl sulfoxide-d ₆	99.96 atom % D
156914-25G	Dimethyl sulfoxide-d ₆	99.96 atom % D
156914-10X0.25ML	Dimethyl sulfoxide-d ₆	99.96 atom % D
417939-5ML	Dimethyl sulfoxide-d ₆	99.96 atom % D, contains 0.03 % (v/v) TMS
417939-25ML	Dimethyl sulfoxide-d ₆	99.96 atom % D, contains 0.03 % (v/v) TMS
417939-10X0.75ML	Dimethyl sulfoxide-d ₆	99.96 atom % D, contains 0.03 % (v/v) TMS
Ethanol		
489336-5G	Ethanol-1,1,2,2,2-d ₅	99.5 atom % D
329347-1G	Ethanol-2,2,2-d ₃	99 atom % D
186414-1G	Ethanol-d ₆	anhydrous, ≥99.5 atom % D
186414-5G	Ethanol-d ₆	anhydrous, ≥99.5 atom % D
151904-100G	Ethanol-OD	≥99.5 atom % D
151904-25G	Ethanol-OD	≥99.5 atom % D
452556-100ML	Ethanol-OD	99 atom % D
452556-25ML	Ethanol-OD	99 atom % D
Ether-d₁₀		
613479-1G	Ether-d ₁₀	99 atom % D
Ethyl acetate-d₈		
522899-1G	Ethyl acetate-d ₈	99.5 atom % D, 99% (CP)
522899-500MG	Ethyl acetate-d ₈	99.5 atom % D, 99% (CP)
522899-5G	Ethyl acetate-d ₈	99.5 atom % D, 99% (CP)
Ethylbenzene-d₁₀		
437344-10ML	Ethylbenzene-d ₁₀	99 atom % D
437344-1ML	Ethylbenzene-d ₁₀	99 atom % D
437344-5ML	Ethylbenzene-d ₁₀	99 atom % D
Ethylene glycol		
343811-1G	Ethylene glycol-(OD) ₂	98 atom % D, 99% (CP)
530549-5G	Ethylene glycol-d ₆	98 atom % D
347442-1G	Ethylene-d ₄ glycol	98 atom % D
Fluorobenzene-d₅		
175803-1G	Fluorobenzene-d ₅	98 atom % D
175803-5G	Fluorobenzene-d ₅	98 atom % D
Heptane-d₁₆		
303011-1G	Heptane-d ₁₆	99 atom % D
Hexafluoroacetone deuterate		
534773-1G	Hexafluoroacetone deuterate	≥99.5 atom % D
1,1,1,3,3,3-Hexafluoro-2-propanol-OD		
411302-25G	1,1,1,3,3,3-Hexafluoro-2-propanol-OD	98 atom % D
411302-5G	1,1,1,3,3,3-Hexafluoro-2-propanol-OD	98 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
Hexane-d₁₄		
303003-5G	Hexane-d ₁₄	99 atom % D
303003-1G	Hexane-d ₁₄	99 atom % D
Imidazole-d₄		
437298-1G	Imidazole-d ₄	98 atom % D
2-Iodopropane-d₇		
377023-5G	₂ -Iodopropane-d ₇	98 atom % D, contains copper as stabilizer
377023-1G	₂ -Iodopropane-d ₇	98 atom % D, contains copper as stabilizer
Methanol		
343854-10G	Methanol-d ₃	99.8 atom % D
343854-5G	Methanol-d ₃	99.8 atom % D
343854-1G	Methanol-d ₃	99.8 atom % D
422878-10ML	Methanol-d ₄	99 atom % D
422878-50ML	Methanol-d ₄	99 atom % D
151947-10G-GL	Methanol-d ₄	≥99.8 atom % D
151947-1G	Methanol-d ₄	≥99.8 atom % D
151947-25G	Methanol-d ₄	≥99.8 atom % D
151947-100G	Methanol-d ₄	≥99.8 atom % D
151947-5G	Methanol-d ₄	≥99.8 atom % D
151947-10G-SB	Methanol-d ₄	≥99.8 atom % D
151947-50G-GL	Methanol-d ₄	≥99.8 atom % D
151947-50G-SB	Methanol-d ₄	≥99.8 atom % D
151947-10G-AMP	Methanol-d ₄	≥99.8 atom % D
151947-1KG	Methanol-d ₄	≥99.8 atom % D
441384-10X0.6ML	Methanol-d ₄	≥99.8 atom % D
441384-10X1ML	Methanol-d ₄	≥99.8 atom % D
441384-10X0.75ML	Methanol-d ₄	≥99.8 atom % D
441384-10X0.5ML	Methanol-d ₄	≥99.8 atom % D
569534-50G	Methanol-d ₄	anhydrous, ≥99.8 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
535435-10G	Methanol-d ₄	≥99.96 atom % D, contains 0.03 % (v/v) TMS
611646-25G	Methanol-d ₄	≥99.8 atom % D, contains 0.05 % (v/v) TMS
611646-10G	Methanol-d ₄	≥99.8 atom % D, contains 0.05 % (v/v) TMS
151939-100G	Methanol-OD	99.5 atom % D
151939-25G	Methanol-OD	99.5 atom % D
550574-25G	Methanol-OD	99 atom % D
550574-100G	Methanol-OD	99 atom % D
Methylcyclohexane-d₁₄		
306053-1G	Methylcyclohexane-d ₁₄	99.5 atom % D
306053-5G	Methylcyclohexane-d ₁₄	99.5 atom % D
Nitrobenzene-d₅		
151955-5G	Nitrobenzene-d ₅	99.5 atom % D
151955-25G	Nitrobenzene-d ₅	99.5 atom % D
151955-10G	Nitrobenzene-d ₅	99.5 atom % D
Nitromethane-d₃		
151963-10X1ML	Nitromethane-d ₃	99 atom % D
151963-25G	Nitromethane-d ₃	99 atom % D
151963-5G	Nitromethane-d ₃	99 atom % D
151963-10G	Nitromethane-d ₃	99 atom % D
269867-5G	Nitromethane-d ₃	≥99 atom % D, contains 1 % (v/v) TMS
Octane-d₁₈		
151971-5G	Octane-d ₁₈	98 atom % D
151971-1G	Octane-d ₁₈	98 atom % D
Pentafluorophenol-OD		
411957-5G	Pentafluorophenol-OD	98 atom % D

Excellent chemical purity and highest isotopic enrichment

Cat. No.	Product Description	Isotopic Purity and Internal Standards
569534-10X1ML	Methanol-d ₄	anhydrous, ≥99.8 atom % D
417653-50G	Methanol-d ₄	≥99.8 atom % D, contains 1 % (v/v) TMS
417653-10G	Methanol-d ₄	≥99.8 atom % D, contains 1 % (v/v) TMS
439029-10G	Methanol-d ₄	≥99.8 atom % D, contains 0.1 % (v/v) TMS
439029-50G	Methanol-d ₄	≥99.8 atom % D, contains 0.1 % (v/v) TMS
444758-10X0.25ML	Methanol-d ₄	99.96 atom % D
444758-10X0.75ML	Methanol-d ₄	99.96 atom % D
444758-1G-AMP	Methanol-d ₄	99.96 atom % D
444758-10X1ML	Methanol-d ₄	99.96 atom % D
444758-10X0.5ML	Methanol-d ₄	99.96 atom % D
444758-10G-AMP	Methanol-d ₄	99.96 atom % D

Cat. No.	Product Description	Isotopic Purity and Internal Standards
Pentane-d₁₂		
490482-5G	Pentane-d ₁₂	98 atom % D
2-Propanol-1,1,1,3,3,3-d₆		
392898-1G	2-Propanol-1,1,1,3,3,3-d ₆	99 atom % D
2-Propanol-d₈		
175897-25G	2-Propanol-d ₈	99.5 atom % D
175897-5G	2-Propanol-d ₈	99.5 atom % D
Pyridine-d₅		
532975-25G	Pyridine-d ₅	≥99.5 atom % D
532975-10G	Pyridine-d ₅	≥99.5 atom % D
532975-10X1ML	Pyridine-d ₅	≥99.5 atom % D
532975-10X0.6ML	Pyridine-d ₅	≥99.5 atom % D
532975-10X0.5ML	Pyridine-d ₅	≥99.5 atom % D
532967-5G	Pyridine-d ₅	≥99.5 atom % D, contains 0.03 % (v/v) TMS
532967-10X1ML	Pyridine-d ₅	≥99.5 atom % D, contains 0.03 % (v/v) TMS

Cat. No.	Product Description	Isotopic Purity and Internal Standards
532967-25G	Pyridine-d ₅	≥99.5 atom % D, contains 0.03 % (v/v) TMS
532967-10X0.6ML	Pyridine-d ₅	≥99.5 atom % D, contains 0.03 % (v/v) TMS
177970-10X0.75ML	Pyridine-d ₅	≥99.96 atom % D
177970-10X0.5ML	Pyridine-d ₅	≥99.96 atom % D
1,1,2,2-Tetrachloroethane-d₂		
358703-1G	1,1,2,2-Tetrachloroethane-d ₂	≥99.5 atom % D
358703-5G	1,1,2,2-Tetrachloroethane-d ₂	≥99.5 atom % D
358703-25G	1,1,2,2-Tetrachloroethane-d ₂	≥99.5 atom % D
425362-5G	1,2-Dibromoethane-d ₄	99 atom % D
425362-25G	1,2-Dibromoethane-d ₄	99 atom % D
331511-1G	1,2-Dichlorobenzene-d ₄	98 atom % D
331511-10G	1,2-Dichlorobenzene-d ₄	98 atom % D
331511-5G	1,2-Dichlorobenzene-d ₄	98 atom % D
396540-5G	1,2-Dichloroethane-d ₄	99 atom % D
396540-1G	1,2-Dichloroethane-d ₄	99 atom % D
186406-1G	1,4-Dioxane-d ₈	≥99 atom % D
302996-1G	1-Butanol-d ₁₀	99 atom % D
Tetrahydrofuran-d₈		
184314-5G	Tetrahydrofuran-d ₈	≥99.5 atom % D
184314-1G	Tetrahydrofuran-d ₈	≥99.5 atom % D
184314-10G	Tetrahydrofuran-d ₈	≥99.5 atom % D
441406-10X1.0ML	Tetrahydrofuran-d ₈	≥99.5 atom % D
441406-10ML	Tetrahydrofuran-d ₈	≥99.5 atom % D
441406-5X.5ML	Tetrahydrofuran-d ₈	≥99.5 atom % D
441406-10X.75ML	Tetrahydrofuran-d ₈	≥99.5 atom % D
269891-10G	Tetrahydrofuran-d ₈	≥99.5 atom % D, contains 1 % (v/v) TMS
269891-1G	Tetrahydrofuran-d ₈	≥99.5 atom % D, contains 1 % (v/v) TMS
437727-1G	Tetrahydrofuran-d ₈	≥99.5 atom % D, contains 0.03 % (v/v) TMS
437727-10G	Tetrahydrofuran-d ₈	≥99.5 atom % D, contains 0.03 % (v/v) TMS
437727-5G	Tetrahydrofuran-d ₈	≥99.5 atom % D, contains 0.03 % (v/v) TMS

Cat. No.	Product Description	Isotopic Purity and Internal Standards
Toluene-d₈		
151998-25G	Toluene-d ₈	99 atom % D
151998-10G	Toluene-d ₈	99 atom % D
151998-50G	Toluene-d ₈	99 atom % D
233382-10X1ML	Toluene-d ₈	99.96 atom % D
233382-10X0.5ML	Toluene-d ₈	99.96 atom % D
233382-1G	Toluene-d ₈	99.96 atom % D
269875-25G	Toluene-d ₈	99 atom % D, contains 1 % (v/v) TMS
434388-10X1ML	Toluene-d ₈	99.6 atom % D
434388-10X0.75ML	Toluene-d ₈	99.6 atom % D
434388-5G	Toluene-d ₈	99.6 atom % D
434388-10G	Toluene-d ₈	99.6 atom % D
434388-25G	Toluene-d ₈	99.6 atom % D
471399-5G	Toluene-d ₈	99 atom % D, contains 0.03 % (v/v) TMS
471399-25G	Toluene-d ₈	99 atom % D, contains 0.03 % (v/v) TMS
570710-50G	Toluene-d ₈	anhydrous, 99.6 atom % D
Trifluoroacetic acid-d		
152005-10X0.5ML	Trifluoroacetic acid-d	99.5 atom % D
152005-10G	Trifluoroacetic acid-d	99.5 atom % D
152005-5G	Trifluoroacetic acid-d	99.5 atom % D
152005-10X1ML	Trifluoroacetic acid-d	99.5 atom % D
152005-25G	Trifluoroacetic acid-d	99.5 atom % D
152005-10X0.75ML	Trifluoroacetic acid-d	99.5 atom % D
2,2,2-Trifluoroethanol		
426237-25ML	2,2,2-Trifluoroethanol-OD	99 atom % D
426237-5ML	2,2,2-Trifluoroethanol-OD	99 atom % D
612197-1G	2,2,2-Trifluoroethanol-1,1-d ₂	99.5 atom % D
396532-1G	2,2,2-Trifluoroethanol-d ₃	≥99.5 atom % D
396532-5G	2,2,2-Trifluoroethanol-d ₃	≥99.5 atom % D
m-Xylene-d₁₀		
175919-5G	m-Xylene-d ₁₀	98 atom % D
o-Xylene-d₁₀		
175900-5G	o-Xylene-d ₁₀	99 atom % D
p-Xylene-d₁₀		
175927-5G	p-Xylene-d ₁₀	99 atom % D

TraceCERT® Certified Reference Materials for quantitative NMR

Quantitative NMR is increasingly used in Pharmaceutical and Chemical Industry as an efficient tool to quantify organic molecules. Most commonly, proton NMR is applied. However, the implementation of qNMR in new fields of application (e.g. metabolomics, biomarker discovery, physiological pathways) brings along more complex molecules and systems, thus making the usage of ^1H -qNMR challenging. The use of other NMR active nuclei, namely ^{31}P or ^{19}F can be a better option in such cases. We offer a comprehensive range of more than 20 internal qNMR standards for ^1H , ^{31}P and ^{19}F nuclei with different chemical shifts and suitable for a variety of different solvents as described on the following pages. To access the products and to get more detailed technical information about the use of quantitative NMR standards please refer to our website SigmaAldrich.com/qnmr where you can also access various publications and articles including our qNMR brochure giving valuable insight and providing tips how to use quantitative NMR.

In 2009, Sigma-Aldrich Switzerland started to apply qNMR under ISO/ IEC 17025 and ISO Guide 34 accreditation to manufacture organic certified reference materials (CRMs). Using a set of internal qNMR standards we built up a considerable portfolio of CRMs for chromatography. So far, more than 200 products are available including pesticides, polyaromatic hydrocarbons (PAH), phenols, plasticizers, cosmetics, antibiotics, air monitoring substances, amino acids, organic pollutants, natural substances and fatty acids. We continuously refined and optimized our qNMR workflows, also taking into account valuable input from CRM users. In the following, we will share some details that will help you to use your NMR instrument for quantitative measurements, achieving maximum accuracy and reproducibility.

Quantification of Organic Compounds

Since most analytical techniques are compound dependent, reliable quantification of organic material is a very challenging task. For example, using HPLC with UV, DAD or fluorescence detection always requires a traceable reference of the very same compound. However, for most organic compounds, no reliable reference material is available. Therefore, the content of an organic material is usually determined by measuring all potential impurities (such as related compounds, water, residual solvents and inorganic impurities) and calculating the content by subtracting the impurity values from a total of 100%.

An alternative to this laborious method is quantitative NMR which, being a relative primary method and having the big advantage, that organic compounds of different chemical structure can be compared quantitatively.

Different referencing techniques have been tested for qNMR, internal as well as external. We usually prefer the internal reference method (**Figure 1**) although the external standard method also certainly has its advantages, e.g., easier recovery of the analyte material, which may be of importance if very expensive material is analyzed. However, with the internal standard method, much higher precision and lower uncertainties can be achieved. Once the materials have been weighed into the same vial, the ratio of analyte and reference stays the same. In contrast to the external standard method, the amount of added solvent, and hence the concentration of the solution, is not critical for the quantification calculation.

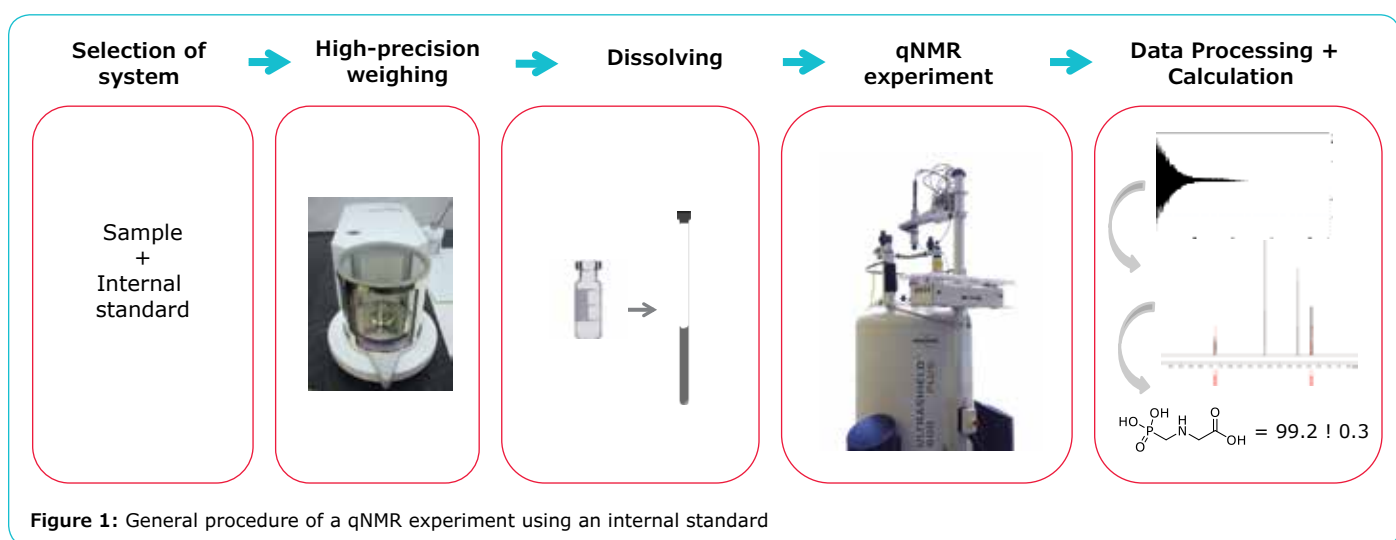


Figure 1: General procedure of a qNMR experiment using an internal standard

Standards for ¹H quantitative NMR, TraceCERT®

PN	Substance	D ₂ O			CDCl ₃			DMSO-d ₆			CD ₃ OD			CD ₃ CN		
		δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)
01380	Ethylene carbonate	4.5	5.5	>250	4.5	7.0	>250	4.5	2.7	>250	4.5	5.3	>250	4.5	2.0	>250
03826	Calcium formate	8.4	3.4	>250	-	-	-	-	-	-	-	-	-	-	-	-
06185	Benzoic acid	-	-	-	8.1 7.7 7.5	3.7 4.0 3.4	150	8.0 7.6 7.5	3.3 3.7 3.0	>250	8.0 7.6 7.5	4.3 4.4 3.9	>250	8.0 7.7 7.5	4.5 2.5 2.6	80
06856	Duroquinone	-	-	-	2.0	3.3	>250	2.0	3.1	30	2.0	4.0	10	2.0	5.4	20
07038	Dimethyl terephthalate	-	-	-	8.1 4.0	3.6 1.8	160	8.1 3.9	2.9 1.1	20	8.1 4.0	4.4 2.4	4.0	8.1 3.9	4.9 2.6	20
14659	Potassium hydrogen phthalate	7.5	2.5	>250*	-	-	-	-	-	-	7.9	2.5	5.0	-	-	-
40384	1,2,4,5-Tetrachloro-3-nitro-benzene	-	-	-	7.8	10.7	>250	8.5	12.6	>250	8.1	6.4	30	8.0	9.6	10
41867	Dimethyl sulfone	3.0	2.9	>250	3.0	2.7	80	3.0	2.4	>250	3.0	3.3	40	2.9	2.6	>250
42582	Ethyl 4 (dimethyl-amino) benzoate	-	-	-	7.9 6.7 4.3 3.1 1.4	4.1 2.4 2.8 2.3 2.6	>250	7.8 6.7 4.2 3.0 1.3	2.7 1.4 1.8 1.5 2.1	>250	7.9 6.7 4.3 3.0 1.4	4.3 2.4 2.9 2.8 2.7	>250	7.9 6.7 4.3 3.0 1.3	5.6 3.7 4.1 3.6 3.6	50
50409	Thymol	-	-	-	7.1 6.8 6.6 3.2 2.3 1.3	3.8 4.5 4.8 4.3 3.1 1.9	>250	7.0 6.7 6.5 3.1 2.2 1.1	2.0 2.2 2.6 2.3 2.0 0.9	>250	7.0 6.6 6.5 3.2 2.2 1.2	3.7 4.4 5.8 4.0 2.8 1.8	>250	7.1 6.7 6.6 3.2 2.2 1.2	4.9 5.7 5.7 5.2 3.5 2.4	>250
55177	Benzyl benzoate	-	-	-	5.4 8.1	4.3 13.6	>250	5.4 8.0	2.2 8.8	>250	5.4 8.1	2.0 3.9	>250	5.4 8.1	1.0 2.4	>250
74599	1,3,5-Trimethoxybenzene	-	-	-	6.1 3.8	4.7 2.2	>250	6.1 3.7	3.2 1.4	>250	6.1 3.8	4.8 2.7	>250	6.7 4.3	3.4 2.6	>250
74658	1,2,4,5-Tetramethyl-benzene	-	-	-	7.0 2.2	6.1 4.0	>250	6.9 2.1	4.7 2.6	10	6.9 2.2	5.9 4.3	10	6.9 2.2	7.7 5.0	50
89151	Dimethylmalonic acid	1.3	1.0	100	-	-	-	1.3	0.7	250	1.4	1.0	250	1.4	2.0	30
92816	Maleic acid	6.3	6.1	>250	-	-	-	6.3	3.0	>250	6.3	3.9	20	6.4	6.2	20
93074	Pentachlorobenzene	-	-	-	7.6	5.8	>250	8.2	13.8	15	7.9	8.2	9	7.8	12.2	15
94681	Methyl 3,5-dinitrobenzoate	-	-	-	9.3 9.2 4.1	8.0 6.1 2.6	>250	9.1 8.9 4.0	8.1 8.1 1.6	100	-	-	-	9.1 9.0 4.0	9.0 8.2 3.5	>250

Standards for ³¹P quantitative NMR, TraceCERT®

PN	Substance	D ₂ O			CDCl ₃			DMSO-d ₆			CD ₃ OD			CD ₃ CN		
		δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)
05498	Triphenyl phosphate	-	-	-	-17.7	2.7	>250	-17.3	1.2	>100	-17.5	3.1	>50	-17.0	4.3	>250
92214	Potassium phosphate monobasic	0.08	8.0	>250	-	-	-	-	-	-	-	-	-	-	-	-
96708	Phosphonoacetic acid	15.7	4.6	>250	-	-	-	14.9	1.5	>250	17.7	2.9	>250	-	-	-

Standards for ¹⁹F quantitative NMR, TraceCERT®

PN	Substance	D ₂ O			CDCl ₃			DMSO-d ₆			CD ₃ OD			CD ₃ CN		
		δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)	δ (ppm)	T1 (s)	s (mg/mL)
07563	4,4'-Difluoro-benzophenone	-	-	-	-105.8	2.4	>250	-106.5	1.4	150	-108.1	2.8	30	-108.3	2.3	140
53396	2,4-Dichloro-benzotrifluoride	-	-	-	-62.5	2.3	>250	-61.2	1.2	>250	-65.4	3.3	>250	-63.0	2.9	>250
80730	2-Chloro-4-fluorotoluene	-	-	-	-115.8	4.4	>250	-115.3	3.3	>250	-117.7	4.8	>250	-117.3	4.7	>250

Table 1: This table gives an overview of spectral shifts, solubility and relaxation times in different solvents of 17 ¹H-qNMR standards, 3 ³¹P-qNMR standards and 3 ¹⁹F-qNMR standards. They cover the entire spectral range, allowing the quantification of almost any organic molecule by qNMR. Solubility tests were done at room-temperature using commercially available NMR solvents. Tests were performed starting from 1 mg/mL up to 250 mg/mL (mg CRM/ml solvent). ₁ relaxation times were recorded for the CRM only (c ≈ 5 to 20 mg/mL at 25 °C), but may vary in the mixture. Therefore it is recommended to check the ₁ delay prior to the qNMR experiment. (*with DCl or NaOD).

Following some basic rules, and with a good balance, it is possible to routinely perform qNMR measurements with 0.5 - 1% measurement uncertainty. Please refer to our website to get more detailed information and tips for a successful qNMR measurement.

In order to minimize the increments by the internal reference, Sigma-Aldrich certified its TraceCERT® CRMs under optimal conditions, so that their uncertainties down to 0.1% would not dominate the overall result, when applying the CRM. Please find a list of the currently available TraceCERT® CRMs for qNMR including some key parameters like typical chemical shifts, relaxation times and solubilities in **Table 1**.

We are currently providing a series of 17 different certified reference materials (CRMs) designed for use in ¹H qNMR experiments and each 3 different CRMs for the use in ¹⁹F or ³¹P qNMR experiments. All products are either traceable to primary material from NIST (National Institute of Standards and Technology) or NMIJ (National Metrology Institute of Japan) and are produced under ISO/IEC 17025 and ISO 17034 (formerly ISO Guide 34) double accreditation.

Figure 2 shows the first page of an example certificate, highlighting the most important features such as traceability statement, expiration date, certified value and associated uncertainty. Since the portfolio is continuously expanding, we also recommend viewing the most current product list on **SigmaAldrich.com/qnmr**.

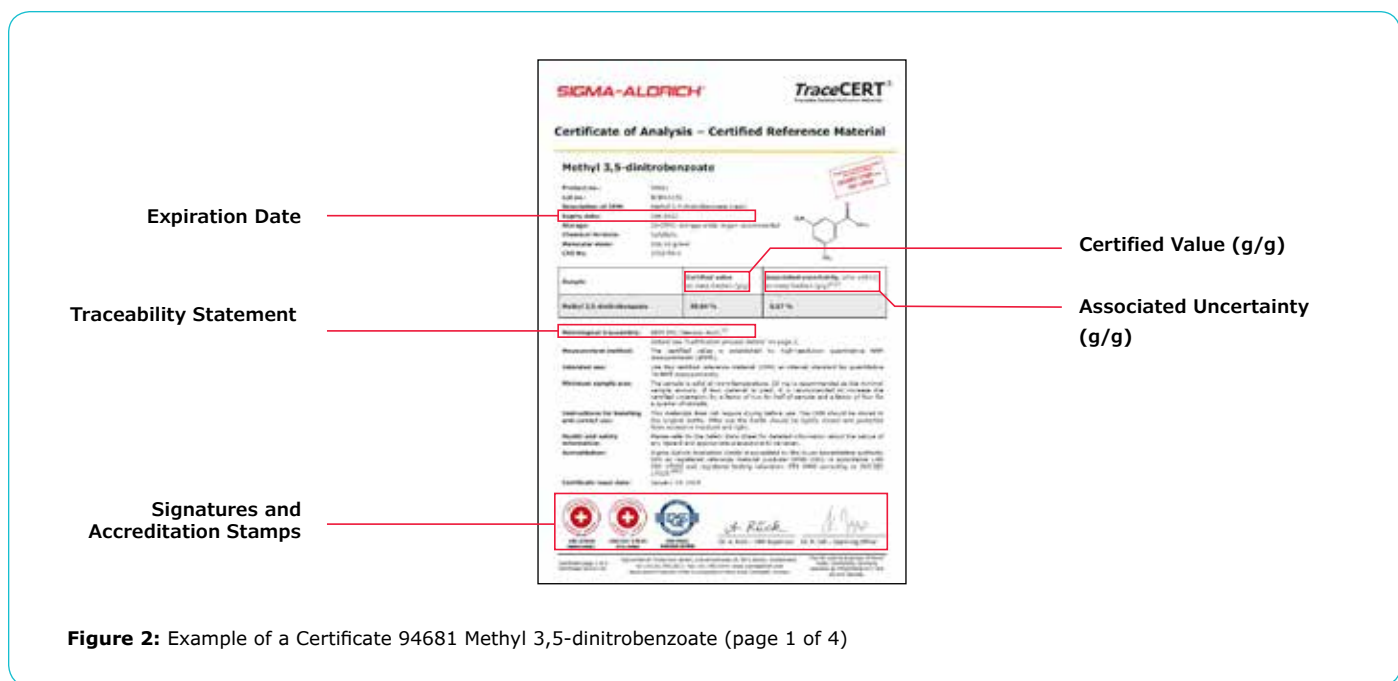


Figure 2: Example of a Certificate 94681 Methyl 3,5-dinitrobenzoate (page 1 of 4)

NMR Reference Standards for Calibration

We also offer a wide range of NMR reference standards needed to verify instrument performance. These standards allow the user to monitor key operational parameters such as pw90, sensitivity, resolution and line shape. We have the combined knowledge and expertise to supply the standards required to ensure proper instrumental qualification and performance

verification. Merck also provides NMR reference standards for Fluorine Sensitivity and Phosphorus Sensitivity verification and other applications. The NMR Reference Standards are available in several tube sizes from 3 mm x 8" to 10 mm x 8". NMR continues to find new and exciting applications to solve today's analytical challenges.

Cat.No.	Concentration	NMR Tube O.D. x L	Application
611778	1% Chloroform in acetone-d6 (99.9 atom % D)	3 mm x 8 in.	1H line shape
487163	1% Chloroform in acetone-d6 (99.9 atom % D)	5 mm x 8 in.	1H line shape
487759	5% Chloroform in acetone-d6 (99.9 atom % D)	5 mm x 8 in.	1H line shape
551341	5% Ethylbenzene, 2% TMS in chloroform-d (99.8 atom % D)	5 mm x 8 in.	1H sensitivity
487112	0.1% Ethylbenzene, 0.01% TMS in chloroform-d (99.8 atom % D)	3 mm x 8 in.	1H sensitivity
487104	0.1% Ethylbenzene, 0.01% TMS in chloroform-d (99.8 atom % D)	5 mm x 8 in.	1H sensitivity
612065	0.1% Ethylbenzene, 0.01% TMS in chloroform-d (99.8 atom % D)	10 mm x 8 in.	1H sensitivity
487147	1% 1,2-Dichlorobenzene in acetone-d6 (99.9 atom % D)	5 mm x 8 in.	1H resolution
551333	12% Tetramethylsilane in chloroform-d (99.8 atom % D)	5 mm x 8 in.	1H reference
487139	0.1 mg/mL Gadolinium(III) chloride, 0.1% DSS, 1.0% water in deuterium oxide (99.9 atom % D)	5 mm x 8 in.	Auto test
612073	40% 1,4-Dioxane, 5 mg/mL chromium(III) acetylacetonate in benzene-d6 (99.6 atom % D)	3 mm x 8 in.	13C PW90, sensitivity
611905	40% 1,4-Dioxane, 5 mg/mL chromium(III) acetylacetonate in benzene-d6 (99.6 atom % D)	5 mm x 8 in.	13C PW90, sensitivity
551368	40% 1,4-Dioxane in benzene-d6 (99.6 atom % D)	5 mm x 8 in.	13C PW90, sensitivity
551384	90% Formamide in dimethyl sulfoxide-d6 (99.9 atom % D)	5 mm x 8 in.	15N sensitivity
487155	1% Iodomethane-13C, 0.2% chromium(III) acetylacetonate, 1% trimethyl phosphite in chloroform-d (99.8 atom % D)	3 mm x 8 in.	Indirect detection test
551406	0.05% α,α,α -Trifl uorotoluene in benzene-d6 (99.6 atom % D)	5 mm x 8 in.	19F sensitivity
611921	25% Hexamethyldisiloxane in benzene-d6 (99.6 atom % D)	5 mm x 8 in.	29Si sensitivity
551392	0.0485M Triphenyl phosphate in chloroform-d (99.8 atom % D)	5 mm x 8 in.	31P sensitivity
611735	0.0485M Triphenyl phosphate in chloroform-d (99.8 atom % D)	10 mm x 8 in.	31P sensitivity
612162	85% Phosphoric acid in deuterium oxide (99.9 atom % D)	5 mm x 8 in. (Coaxial)	31P sensitivity

Stable Isotope Labeled Reagents and Building Blocks

Merck is the market leader in organic and inorganic stable isotope labeled chemicals, providing an innovative and comprehensive range of building blocks, reagents and advanced materials for chemical synthesis, medicinal chemistry and materials science. Look to us for high-quality reagents and intermediates for your research product needs.

ISOTEC® Stable Isotopes specializes in the custom synthesis of labeled compounds. employs a highly trained group of stable isotope scientists. The stable

isotope R&D group routinely engages in complex, multi-step syntheses of a wide range of complex compounds with a variety of labeling patterns.

Our ISO 9001 facilities give us the capability to perform synthesis ranging from milligram to kilogram quantities. In addition, we are able to offer a variety of product grades from standard research reagents through full GMP API manufacturing.

Cat. No.	Product	Isotopic Purity
Acids and Bases		
282022	Acetic acid- ¹³ C ₂	99 atom % ¹³ C
283835	Bromoacetic acid- ¹³ C ₂	99 atom % ¹³ C
543047	Deuterium chloride solution, 35 wt. % in D ₂ O	99 atom % D
279404	Formic acid- ¹³ C, 95 wt. % in H ₂ O	99 atom % ¹³ C
426229	Formic acid-d ₂ , 95 wt. % in D ₂ O	98 atom % D
609323	Nitric- ¹⁵ N acid solution, ~10 N in H ₂ O	98 atom % ¹⁵ N
Inorganic Salts		
366536	Ammonium nitrate- ¹⁵ N	98 atom % ¹⁵ N
297038	Potassium cyanide- ¹³ C	99 atom % ¹³ C
490539	Potassium cyanide- ¹³ C, ¹⁵ N	99 atom % ¹³ C, 98 atom % ¹⁵ N
190020	Sodium cyanoborodeuteride	96 atom % D
279412	Sodium formate- ¹³ C	99 atom % ¹³ C
364606	Sodium nitrate- ¹⁵ N	98 atom % ¹⁵ N
490814	Sodium nitrite- ¹⁵ N	98 atom % ¹⁵ N
Aldehydes		
531227	Acetaldehyde- ¹³ C ₂	99 atom % ¹³ C
489417	Formaldehyde- ¹³ C solution, 20 wt. % in H ₂ O	99 atom % ¹³ C
596388	Formaldehyde- ¹³ C, d ₂ solution, 20 wt. % in D ₂ O	98 atom % D, 99 atom % ¹³ C
492620	Formaldehyde-d ₂ solution, ~20 wt. % in D ₂ O	99 atom % D
485233	4-Hydroxybenzaldehyde-1- ¹³ C	99 atom % ¹³ C
Buffers		
473111	Lithium- ⁶ Li ₂ carbonate	95 atom % ⁶ Li
372382	Sodium bicarbonate- ¹³ C	98 atom % ¹³ C
490768	Sodium carbonate- ¹³ C	99 atom % ¹³ C
486248	Tris-d ₁₁ solution, 1 M in D ₂ O	98 atom % D

Cat. No.	Product	Isotopic Purity
Alcohols		
451835	Bisphenol A-d ₁₆	98 atom % D
440671	1,1,1,3,3,3-Hexafluoro-2-propanol-d ₂	99 atom % D
490504	Phenol- ¹³ C ₆	99 atom % ¹³ C
Amines		
485497	Aniline- ¹³ C ₆	99 atom % ¹³ C
617490	Benzidine-(rings-d ₈)	98 atom % D
488208	Benzylamine- ¹⁵ N	98 atom % ¹⁵ N
488828	Diethyl-d ₁₀ -amine	98 atom % D
676152	Formamidine- ¹³ C, d, ¹⁵ N ₂ acetate	98 atom % D, 98 atom % ¹⁵ N, 99 atom % ¹³ C
492787	Hydrazine- ¹⁵ N ₂ monohydrate	98 atom % ¹⁵ N
296007	Methyl- ¹³ C, d ₃ -amine hydrochloride	99 atom % D, 99 atom % ¹³ C
176001	Methyl-d ₃ -amine hydrochloride	99 atom % D
793752	p-Toluidine-(phenyl- ¹³ C ₆)	99 atom % ¹³ C
486353	tert-Butyl-d ₉ -amine	98 atom % D
Halides		
488232	Bromobenzene- ¹³ C ₆	99 atom % ¹³ C
701734	Bromobenzene-d ₅	99 atom % D
293199	Ethyl bromoacetate-1- ¹³ C	99 atom % ¹³ C
283800	Ethyl bromoacetate- ¹³ C ₂	99 atom % ¹³ C
293172	Ethyl bromoacetate-2- ¹³ C	99 atom % ¹³ C
721328	Iodoacetamide- ¹³ C ₂ , 2-d ₂	98 atom % D, 99 atom % ¹³ C
277185	Iodomethane- ¹³ C, contains copper as stabilizer	99 atom % ¹³ C
176036	Iodomethane-d ₃ , contains copper as stabilizer	99.5 atom % D

NMR Tubes and Accessories

Merck is pleased to carry a broad range of NMR tubes and accessories.

- Wilmad®, Norell®, Shigemi and Aldrich NMR tubes
- Septa seals
- Tube caps
- Tube carriers
- Tube labels
- Tube racks

Premium NMR Tubes

These tubes have the tightest specifications for concentricity and camber, ensuring a well-balanced spin rate and the highest shimming quality.

ASTM type 1, Class A thin wall (0.38 mm) borosilicate glass

- **Aldrich® Precision Grade ColorSpec® NMR tubes**
RED color-coded banding and caps, and a MHz rating that allows for instant visual identification.
- **Norell® Select Series™ NMR tubes**
- **Shigemi NMR Microtubes**
Magnetic susceptibility matched to CDCl₃, CD₃OD, DMSO-d₆, or D₂O for best resolution.
- **Wilma d® Precision Grade NMR tubes**

For a complete listing of all available High-Resolution NMR solvents and tubes, visit SigmaAldrich.com/highresnmr

Economy NMR Tubes

NMR tubes for routine applications allow for slightly greater variance of concentricity and camber than precision tubes and have a greater impact on shimming quality. With proper handling and sample preparation, they provide sufficient quality for most standard analysis.

ASTM Type 1, Class B medium wall (0.43 mm) borosilicate glass

- **Aldrich® Thrift Grade ColorSpec® NMR tubes**
YELLOW color-coded banding and caps, and a MHz rating that allows for instant visual identification.
- **Norell® Standard Series NMR tubes**
- **Wilma d® Economy Grade NMR tubes**

For a complete listing of all available Routine NMR solvents and tubes, visit SigmaAldrich.com/routinenmr

High-Throughput NMR Tubes

Type 1, Class A or Class B borosilicate glass. Sold in larger packs (50 or 100) and are disposable.

- **Aldrich® Disposable ColorSpec® NMR tubes**
BLUE color-coded banding and caps, and a MHz rating that allows for instant visual identification.
- **Norell® Standard Series NMR tubes**
- **Wilma d® Economy Grade NMR tubes**
- **Bruker® Match tubes**
- **Norell® Sample Vault Series NMR tubes**

For Bruker® SampleJet™ robotic high-throughput automated systems.

For a complete listing of all available High-Throughput NMR solvents and tubes, visit SigmaAldrich.com/htpnmr

Specialty NMR Tubes and Accessories

Offering a full range of products designed to eliminate the effects of exposure to air and moisture.

- **Wilma d® screw-cap NMR tubes (Z271950)**
Open-top, screw-thread cap with a PTFE/silicone septum that permits addition of reactants or samples via syringe.
- **Wilma d® quick pressure valve NMR tubes (Z562866)**
Designed to simplify the process when performing studies in catalysis, gas-liquid phase reactions, air-sensitive samples and elevated temperature experiments.
- **Wilma d® low pressure/vacuum valved NMR tubes (Z568139)**
A convenient flame-free sealing solution for air-sensitive or volatile liquid samples.
- **Aldrich® 5 mm NMR tube filling manifold (Figures 1 – 3)**
Convenient, easy system for flame sealing air- and moisture-sensitive NMR tubes. Manifolds available in 1 port, 2 port and 5 port configurations.
- **Aldrich® NMR-tube Dewar flask**
Designed for freezing NMR samples in vacuum-thaw work.
- **eVol® electronic syringe NMR starter kit**
The eVol® NMR Edition kit allows accurate manipulation of small volumes of NMR samples enabling dilution directly in the analysis tube — even in tubes with the smallest internal diameters.
- **Precision Seal® rubber septa caps for 5 mm tubes**
Provides a penetration point for cannulation and a dual seal inside the tube.

For a complete listing of all available Air-Sensitive NMR solvents, tubes and NMR accessories, visit SigmaAldrich.com/airsensnmr

TRUST OUR NMR EXPERTISE

For Products that meet your specs

- Deuterated Solvents and Reagents
- Isotec® Stable Isotopes
- Reference Standards
- NMR Tubes and Accessories

[SigmaAldrich.com/nmr](https://www.sigmaaldrich.com/nmr)



Sigma-Aldrich®

Lab & Production Materials

Merck KGaA
Frankfurter Strasse 250
64293 Darmstadt, Germany

SigmaAldrich.com

To place an order or receive technical assistance in Europe, please visit **SigmaAldrich.com/order** or call Customer Service:
Austria: +49 7329 972499 France: 0800 211 408
Germany: +49 7329 972499 Italy: +39 02 3341 7310
Portugal: +351 219 242 555 Spain: +34 916 61 99 77 Option 1
Switzerland: 0848 645 645 United Kingdom: 0870 900 4645

For other countries across Europe, please call:
+44 (0) 115 943 0840

Or visit: **SigmaAldrich.com/offices**
For Technical Service visit: **SigmaAldrich.com/techservice**

© 2019 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. Merck, the vibrant M, Sigma-Aldrich, Aldrich TraceCERT ColorSpec Precision Seal and ISOTEC are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources

2017 - 08821 MK_BR1420EN
01/2019

